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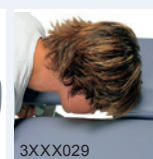
933

851


ACCESSORY: SAUNDERS CLINICAL TRACTION DEVICE

CERVICAL TRACTION WITH THE SAUNDERS CLINICAL TRACTION DEVICE

Cervical traction can be very effective in the treatment of cervical hernias that require widening of the neural foramen or stretching of the soft tissue in the cervical structure. The Saunders clinical traction device makes cervical traction simple and comfortable.

With the Saunders clinical traction device, the traction forces are directed towards the occiput. The forces are not applied to the chin, so there is no pressure on the temporomandibular joint. The Saunders clinical traction device can be adjusted to individual needs so each patient receives the maximum benefit.

The angle of the cervical traction is adjustable with the Saunders clinical traction device. This permits optimal traction as described in the literature.

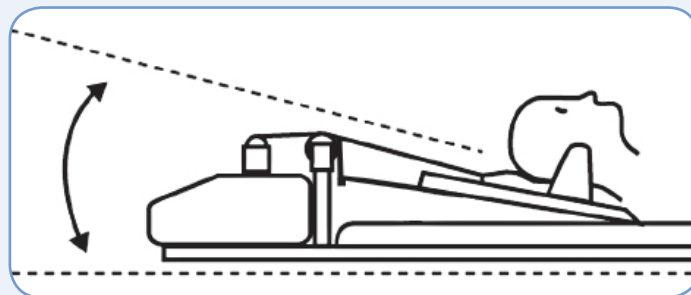The combination of the Eltrac 471 with an Enraf-Nonius traction table and the Saunders clinical traction device is the ideal solution for anyone who is looking for responsible, safe and comfortable cervical traction.

SAUNDERS CLINICAL TRACTION DEVICE

1471808 Saunders clinical traction device including wall bar bracket for the Eltrac 471

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Whether it concerns the treatment of a hernia, degenerative discopathy, or a pinched nerve, spinal decompression provides a solution for prolonged symptoms of pain. Spinal decompression can be achieved in a variety of ways, but traction is one of the easiest and most effective methods. Traction therapy influences the pressure between the vertebrae. Because of the pressure decrease during traction, a mechanical effect on the nerves is generated and a biophysical effect on the circulation and oxygen supply to the local structures. This can result in immediate pain relief, but also a structural improvement in the long term.

The Eltrac 471 is your perfect assistant for traction therapy. The trusted Eltrac 471 has helped many patients for decades on end. This Eltrac 471 has now been improved. The Eltrac 471 is accurate, easy-to-use, real-time and safe!

ACCURATE

The new Eltrac is precise. The forces are exercised on the patient with ultimate precision. The tractive effort in the cord of the Eltrac 471 is constantly measured during traction. The control of the Eltrac compares the measured value to the set value. Any deviations are immediately corrected. The movements of the patient are monitored, such that these have no influence on the actual tractive effort.

The Eltrac is also extremely accurate in reproducing the values of the traction therapy. The actual values are shown, meaning the actual force exercised, rather than indications of a percentage of the motor capacity or a predicted force! And that accuracy is reproduced in the smallest conceivable steps. By zooming in on the display, you can see what is actually happening!

EASY-TO-USE

The new Eltrac 471 is extremely easy-to-use thanks to the touchscreen and the intuitive interface. A clear and well-organised menu will assist you in making therapeutic choices. All essential information will be displayed. Not only while adjusting the settings for the therapy, but also during the therapy itself.

Moreover, the treatment protocols can be saved as well. Every possible combination can be added to the protocol library, even protocols that combine static and intermittent traction in a single session, which was never possible until recently! Static traction alone is often too painful for the patient, particularly in the acute phase when traction can prove to be so effective. In such cases, combining static and intermittent traction in a single session may be the ideal solution. And the new Eltrac 471 offers various forms to choose from in that respect, such as pro- and regressive, cyclic patterns or combinations of these.

REAL-TIME

What you see on the clear and informative display is exactly what is taking place. No averaging later on, but rather real-time. This allows you to monitor the forces actually exercised in the course of the therapy, not a summarised representation, but rather the actual course. This will also provide better insight into the course of the treatment session as a whole. You will be on top of what you need to know thanks to the real-time information.

SAFE

The Eltrac 471 by Enraf-Nonius is the safe and responsible choice. Not only because Enraf-Nonius has been supplying physiotherapists with equipment as early on as 1925, but also because we make up part of the treatment equipment in terms of technological advancement. The traditional equipment for traction is mainly characterised by controlled traction forces. In addition to the controlled contraction force, the new Eltrac 471 also offers clear real-time feedback, the option to save protocols and the user-friendliness that modern technology renders possible. In addition to the option to save your individual protocols, the Eltrac is of course also equipped with pre-programmed clinical protocols (evidence based!). This makes it extremely easy to adjust the settings of the equipment.

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ELTRAC 471

ERGONOMIC DESIGN

A design that focusses on functionality and user-friendliness. Offering the therapist comprehensive insight and a clear interface. Offering the patient the most comfortable fixation method and a relaxed starting position.

THE NEW GENERATION

The Eltrac by Enraf-Nonius is the newest generation of equipment for traction therapy. Thanks to the latest technological possibilities, the most advanced scientific insights and the many years of experience of Enraf-Nonius as developer of therapeutic equipment, you will be guaranteed of having the best equipment with the new Eltrac 471.

TRACTION BENCH CAN ALSO BE USED AS TREATMENT BENCH

The rolling blade components for traction can be easily secured. Your traction bench (type 916 and 926) converts into an ideal treatment bench in this way, adjustable in height and with a folding leg component and head component.

EASY-TO-USE

The new Eltrac is extremely easy-to-use thanks to the touch-screen and the intuitive interface. A clear and well-organised menu will assist you in making therapeutic choices. All essential information will be displayed. The graphic support with the zoom option and direct feedback is unsurpassed.

PRECISE

The forces can be adjusted with the utmost precision. Not a percentage of the motor force, but in precise small steps of 1 N. or 10 N.

PROGRAMMED FOR SPEED

The menu structure is optimized to speed. In most cases you can start a pathology-specific treatment in just 3 steps! Easy, quick and simple.


ROLLING BLADE COMPONENTS

There is maximum comfort for the patient during the traction thanks to the rolling blade components. The forces are exercised fully on the lumbar vertebrae and do not disappear in unpleasant friction between the patient and the bench.

ACCURATE

Feedback of the measured forces exercised and immediate correction of deviation where necessary. In addition, the device precisely adjusts the traction force as quickly as you want.

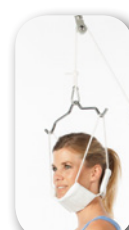MADE BY ENRAF-NONIUS

Made by Enraf-Nonius and, therefore, safe, reliable and durable. A device made by Enraf-Nonius, THE supplier for the physiotherapist worldwide. Where quality, reliability and durability come as standard.


TRACTION BENCH WITH EASY-TO-USE SECURING POINTS

With fixation bars on both the sides as well as the head, the Manumed Traction offers a fixation possibility for every traction. Using the traction belt is the optimal method, but the support of various Velcro straps when applying traction becomes an easy task on a Manumed Traction.


CRYSTAL CLEAR TFT COLOUR SCREEN

A crystal clear, full colour and scratch resistant screen that allows full visibility on all parameters from all viewing angles.

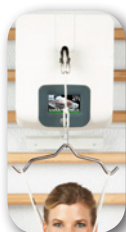PROGRAMMABLE

In addition to the standard treatment protocols, you can also save your personal protocols (including sequential applications) based on your own guidelines. Every conceivable combination - static, pro- and regressive, cyclic or intermittent can be stored as a personal favourite in this way.

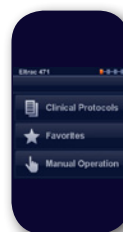TOUCH SCREEN USER INTERFACE

Easy to use by means of a clear menu via a touch screen interface. Extremely intuitive and as easy as a state-of-the-art smart phone.

FLEXI-BENCH FOR AN IDEAL STARTING POSITION DURING LUMBAR TRACTION

The high-low adjustable flexi-bench render it possible to exercise traction with an ideal starting position in which the hip-angle can be bent to 45 to 90 degrees.

AREAS OF APPLICATION OF TRACTION THERAPY

PAIN REDUCTION


REDUCTION OF MUSCLE TONE


INCREASE OF THE RANGE OF MOTION AND FLEXIBILITY

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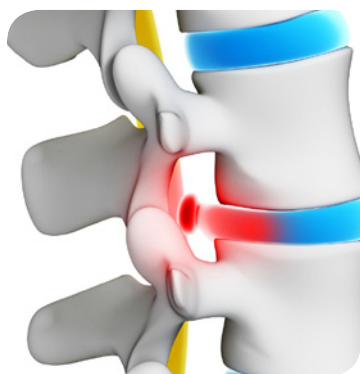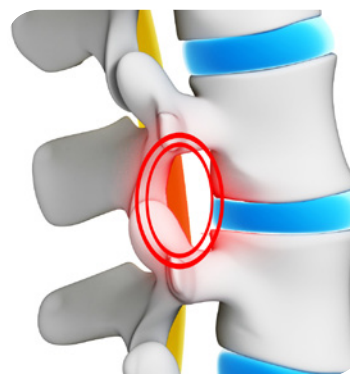
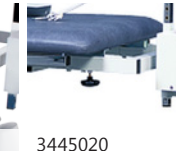
5500926
5500916

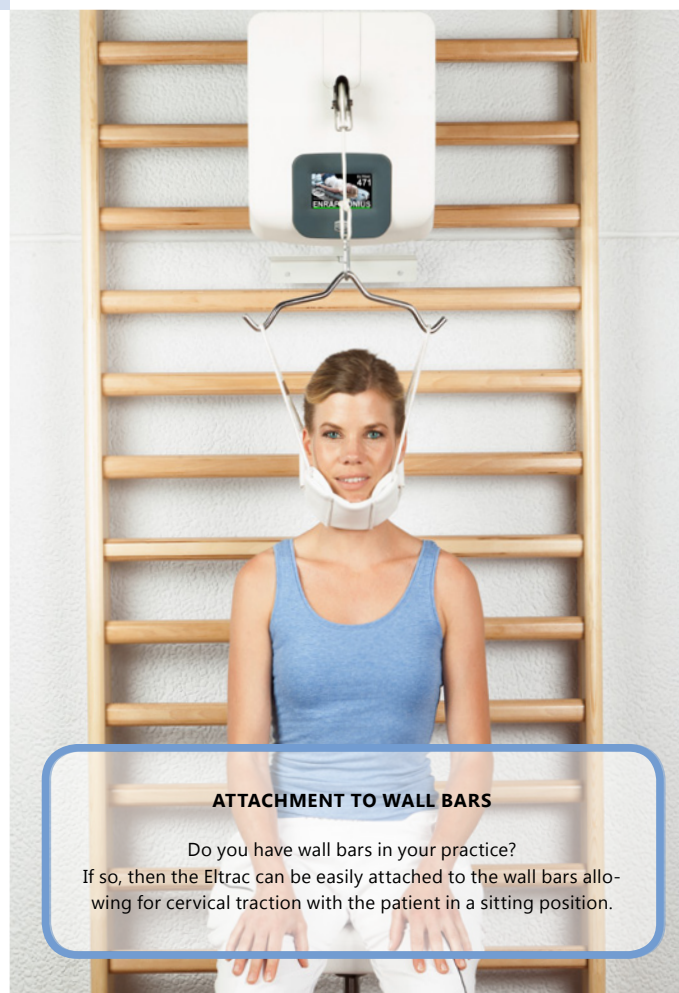ATTACHMENT TO WALL BARS

Do you have wall bars in your practice?
 If so, then the Eltrac can be easily attached to the wall bars allowing for cervical traction with the patient in a sitting position.